



Bellco Glass, Inc.

Biotechnology Manufacturing

Bellco Glass, Inc. presents

Digital R/A Base & 6 Decks, 115V


SKU : 7631-75016

PRODUCT QUALITY

Bellco Glass, Inc. guarantees quality in each and every product that goes out the door, from start to finish. With our strict QC policies and procedures, you can safely bet that you will get the best products Bellco Glass, Inc. has to offer with every single order.

OUR VISION

We will continue to move forward into the future with quality manufacturing, quality raw materials, and constantly redefining and upgrading our internal processes, thus, making us one of the best biotechnology manufacturing companies in the world.

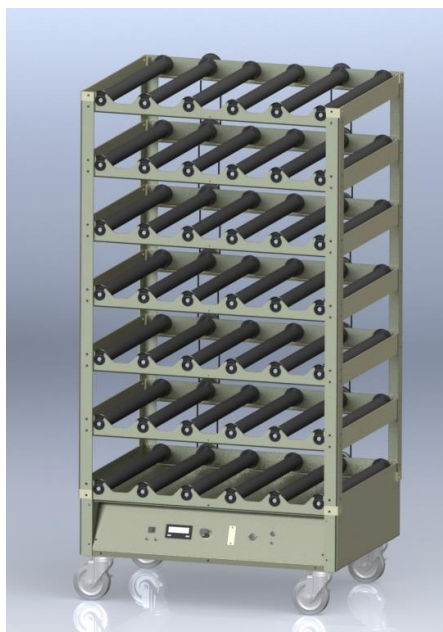

- Bellco Roller Apparatuses are ideal for applications such as the production of mammalian cells, viruses, and cell products
- Modular units are available from 5 to 45 bottle positions and allow easy expansion with modular add on decks
- Rollers have a key design to eliminate slipping and are easily removeable for servicing of bearings, easy remove tool included

Features a solid state control with heavy duty DC motor and a digital LCD display

The use of a brushless motor with controller ensures smooth operation and no maintenance

Front panel controls include a locking speed control knob, and a rotation switch to change from CW/CCW operation

Bellco Glass, Inc. services available for products and customers:

- Technical Support
- Customer Service
- Application Support
- Hardware Support
- Sterilization
- Amber coloring
- Plastic coating



Bellco Glass, Inc.
340 Edrudo Rd.
Vineland, NJ 08360
1-800-257-7043
www.bellcoglass.com



Bellco Glass, Inc.

Biotechnology Manufacturing

Specifications

- 115 +/- 10VAC @ 60 Hz line voltage
- 3 Amp @ 115V, 60 Hz power consumption
- Dimensions: 30"w x 24 3/4"d x 57"h (76 x 63 x 145 cm)
- Number of positions: 35
- Bottle Speed: 0.1 to 4.0 RPM
- Roller Construction: Non-transferring chemical resistant rubber rollers with key shape
- Clearance between Decks for Bottle Placement: 134mm (5 1/4")
- Usable Roller Length: 51 cm (20")
- Construction Material: Aluminum, chemically treated enamel finish
- Bottle Size: Accepts Bellco Glass or Plastic Roller Bottles 110mm to 120mm OD up to 510mm in length
 - 7630-50555 All Position Drive is recommended for use with plastic bottles

TECHNICAL SUPPORT


Technical support is available to all of our customers from inquiry to purchase and through the use of the product. Our service department is dedicated to ensuring the best white glove service to you as possible and to providing quality inspection and review of your concerns, should they arise.

CUSTOMIZATION

A lot of items that come out of Bellco Glass, Inc. go through a customization process. A customer will come to us with a specific need. Our engineering department will review the item and its capabilities to do whatever we can to create a custom item to perform exactly how you need it to. Contact us today for your customized quote.

QUALITY GLASS

Our manufacturing facility uses Type 1 Class A borosilicate glass. This is a type of glass that is of superior thermal and chemical resistance and specifically designed for volumetric glassware, providing you with products that give you longevity of use.

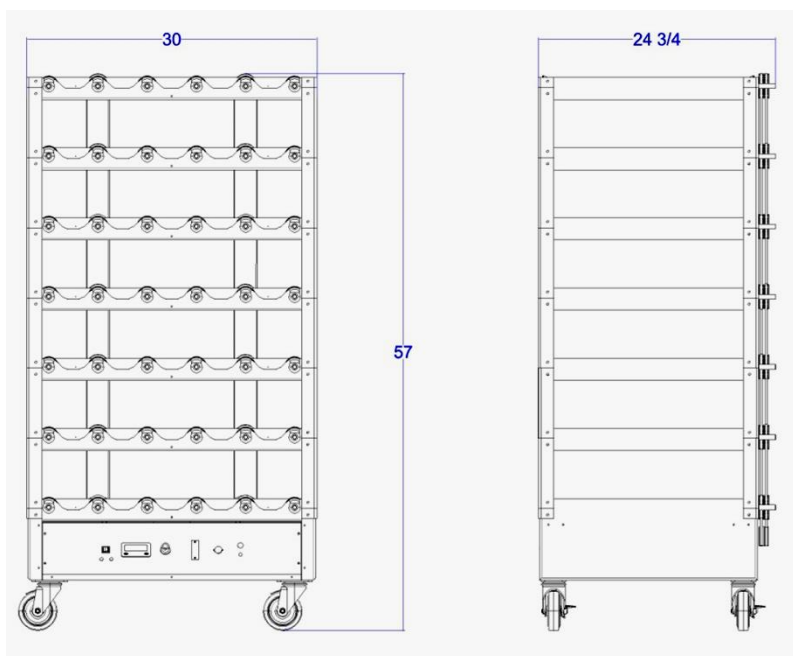

**DISCLAIMER :The information contained in this specification sheet is correct to the best of our knowledge. The recommendations and suggestions therein are made without guarantee or representation as to the results. We recommend that adequate tests be made in your laboratory or plant to determine if this product meets all of your requirements.*



Bellco Glass, Inc.
340 Edrudo Rd.
Vineland, NJ 08360
1-800-257-7043
www.bellcoglass.com