



- The Patented VersaCap® GL45 is designed to withstand autoclaving at 15 PSI and 121°C (250°F) for up to 20 minutes.
- The patented VersaCap® GL45 accommodate a variety of tube adaptors, for venting and fluid transfer.
- The cap is equipped with a stepped barb connection, ensuring compatibility with various tubing sizes and providing a secure, leak-free fit.

Caution

- Over-tightening of caps or connectors can damage threads or compromise the seal.
- Do not exceed the vacuum or pressure limits specified by the manufacturer to prevent damage or accidents.
- Avoid exposing vent filters to liquid or debris, which could block air exchange and compromise performance.
- Change vent filters regularly to maintain sterility and efficiency.
- Ensure all ports and connectors are securely sealed to prevent spills or leaks during use.

The Vactrap™ G Assemblies are designed as modular, flexible, and interchangeable components to upgrade your existing Vactrap™ G system.

The Vactrap™ G Cap Assemblies fit the Standard GL45 Glass Bottle Neck, they come complete with TPE Tubing, PTFE Vent Filter and Plug.

Whether you're setting up a new vacuum trap system or replacing parts in an existing setup, the modular nature of these assemblies simplifies installation while maintaining functionality.

❶ Patented VersaCap® GL45:

Open-top caps made from white polypropylene (PP), providing secure sealing and chemical resistance.


❷ **Foxx Adapter Inserts:** Molded adapters designed specifically for use with glass bottles, ensuring compatibility and ease of use.

❸ **Vent Filter:** 0.2 µm PTFE filter for efficient air exchange, preventing contamination.

❹ **Tubing:** Flexible and durable TPE tubing designed superior performance under vacuum pressure.

❺ **Pinch Clamp:** Durable and adjustable clamp for precise control over tubing flow, ensuring accuracy.



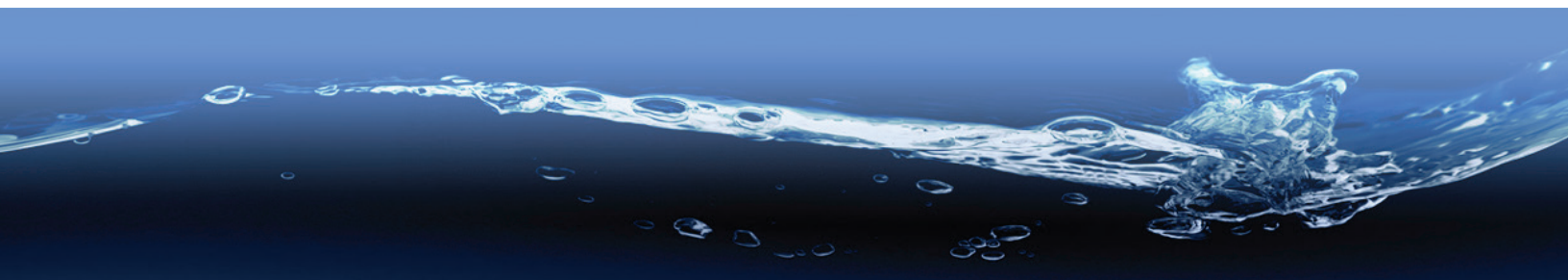

Item	Description	Thread Type	Quantity
300-0037-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Dual Bottle Setup Assembly	GL45	1EA
300-0038-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Single Bottle Setup Assembly	GL45	1EA
300-0039-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Single Bottle Setup Assembly	GL45	1EA
300-0040-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Dual Bottle Setup Assembly	GL45	1EA
300-0041-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Dual Bottle Setup Assembly	GL45	1EA

COMPATIBILITY CHART

Item	Description	Compatible Vacrap System SKU
300-0037-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Dual Bottle Setup Assembly	30R-1101-FLS
		30R-2101-FLS
		30R-2201-FLS
		30R-3101-FLS
300-0038-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Single Bottle Setup Assembly	30R-2001-FLS
300-0039-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Single Bottle Setup Assembly	30R-5001-FLS
300-0040-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Dual Bottle Setup Assembly	30R-5201-FLS
300-0041-FLS	Vacrap™ G Cap Assembly, GL45 Versacap®, PP, TPE Tubing, For Dual Bottle Setup Assembly	302-1001-FLS
		30R-C501-FLS